

Advanced design for complex cases.

Long lasting rings – flash-free results
Drawn-wire nickel titanium produces lasting separating pressure plus ultra-adaptive Soft-Face™ tips for reduced flash.

Enhanced cervical seal
The split-tip Quad wedges combined with the driver tip on the Quad rings produce an independent cervical seal on both sides of back-to-back restorations.

Quadrant workflow
Asymmetrical tip design improves ring fit, matrix stabilization and tooth separation all of which are key for achieving optimal results on challenging multiple-tooth and back-to-back Class II preparations.



Techniques for success
See Quad cases, technique videos, and frequently asked questions by following this link.

purchase options



Quad Sectional Matrix System Intro Kit
QR-KHF-10

- 2x Quad rings (1x each Blue orientation)
- 50x Firm bands (10x each)
- 50x Dead-soft bands (10x each)
- 80x Quad wedges (20x each)
- 10x Varistrip
- 12x Rally polishers
- 1x FXP01 ring forceps

Introductory Price
\$650

Quad Sectional Matrix System Professional Kit
QR-KHF-11

- 4x Quad rings (1x each Blue and Orange orientation)
- 50x Firm bands (10x each)
- 50x Dead-soft bands (10x each)
- 80x Quad wedges (20x each)
- 10x Varistrip
- 12x Rally polishers
- 1x Ring forceps
- 1x Band forceps
- 1x Multi-function instrument

Introductory Price
\$1212



THE LEADER IN MATRIX SYSTEMS

order toll free 888.437.0032 | www.garrisdental.com
Offer valid through 3/31/2025. Cannot be combined with other promotions or discounts unless otherwise specified. Excludes customers with pre-negotiated pricing, including DSOs and buying groups.

PL225QUAD

© 2025 Garrison Dental Solutions, LLC



Quad
Matrix System™

Advanced design
for complex cases


Garrison

Garrison

FROM THE LEADER IN MATRIX SYSTEMS

Go beyond basic Class IIs with Quad.

While it's great for everyday standard Class IIs, the Quad Matrix System is named for its unique ability to streamline more complex cases. Faced with an entire quadrant of Class IIs like this, you need unique material designs and proven techniques to deliver efficient results you can be proud of.

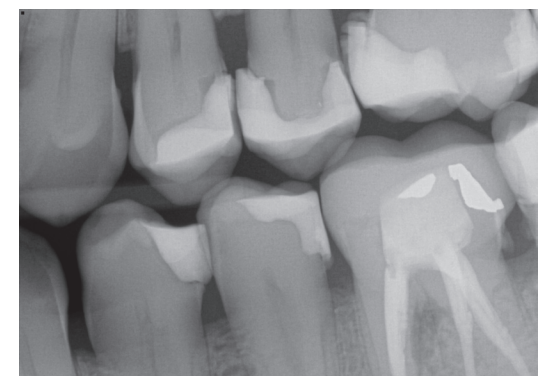

When a case like this sits down in your chair, you need Quad.

ClinicianDesigned–Garrison® Engineered.



When Dr. Matt Burton, DDS, a chair-side dentist, envisioned a better way to tackle complex Class II restorations—like back-to-back cases, deep lesions, and multi-tooth restorations—Garrison's expertise in matrix systems and highly engineered products made it a reality. This innovative system, backed by thousands of hours of clinical field validation, overcomes challenges with precision and efficiency.

This collaboration between clinical vision and engineering mastery will streamline your workflow for all Class IIs.

Quad Matrix System™

A fully integrated sectional matrix system.

Confidently address challenging restorations including:

Back-to-Back Restorations: Simultaneously restore both sides of adjacent Class IIs with confidence. Two matrix bands, one Quad Wedge, and one Quad Ring make for a streamlined workflow on this common procedure.

Deeper Caries: Use a subgingival Firm Band along with a Quad Wedge for a deeper more secure seal enhanced by the Quad Ring driver tip.

Complex Multi-Tooth Restorations: Use combinations of Universal (blue) and Tall (orange) Quad rings to produce ideal contours and contacts on the most challenging cases. We'll show you how!



Advanced design for complex cases.



All-new asymmetrical **Quad Ring** design with specialized 'driver-tip' maximizes Quad Wedge adaptation and cervical seal.




- Adapts the split-tip Quad wedge to both sides of the embrasure simultaneously
- Seals both cervical margins
- Drives deeper for improved matrix band adaptation and tooth separation
- Soft-Face™ design for superior flash reduction!




Maximum adaptation
The right and left driver tips of the Quad system allow you to select the ideal ring for each quadrant and for each restorative challenge.




Tall rings (orange) with both left and right-side driver tips are available for simplified stacking.

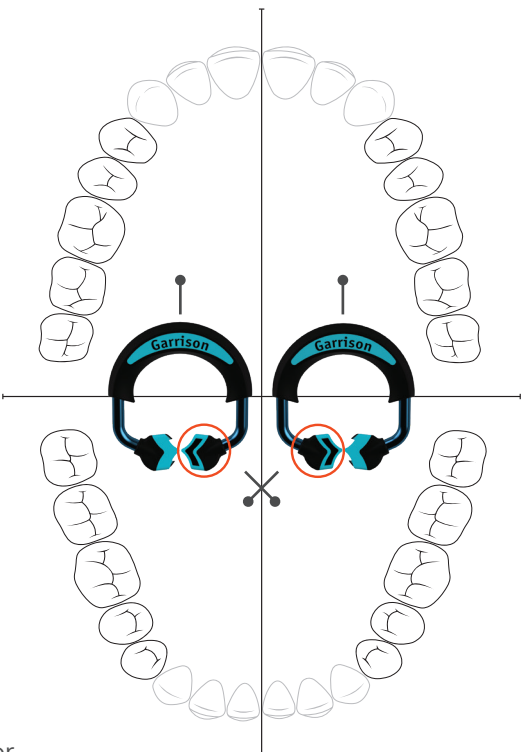

The **Quad Wedge** utilizes a revolutionary split-tip design.

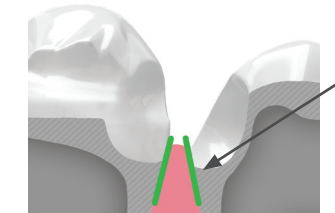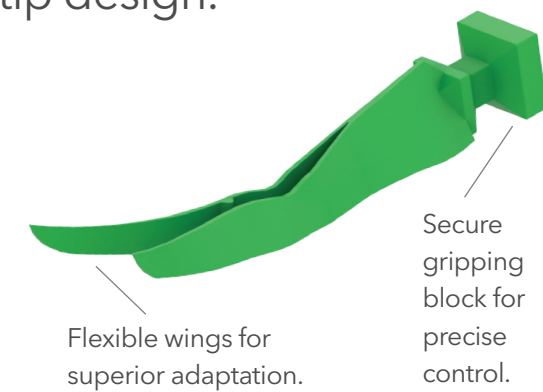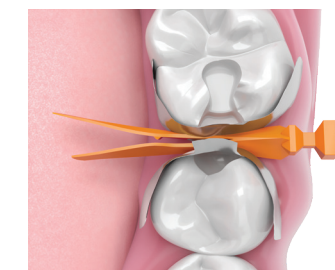


Split tip easily adapts bands to both cervical margins!
Gently slides alongside papilla and seals deeper for a more natural emergence profile— even with mismatched depths.

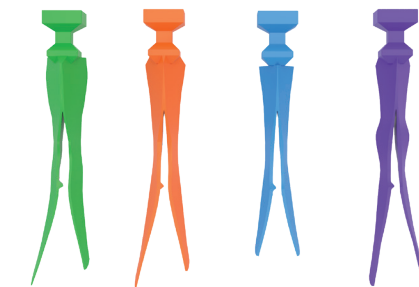


Three sizes plus the special, anatomically-shaped purple!

Seal those concavities!
Quad Wedges include the highly specialized anatomical purple wedge designed for an improved seal when presented with a root concavity.

- Prevents flash and contamination.
- Improves matrix adaptation.

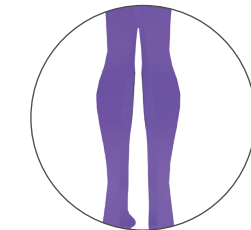


A comprehensive variety of **Matrix Bands** to accommodate a wide range of clinical scenarios.

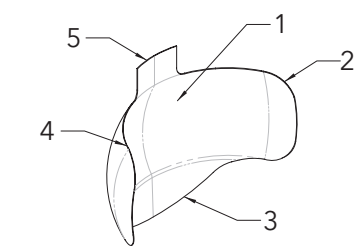

1. Ideal interproximal anatomy built in
2. Ultra thin!
a. Firm Bands 0.0015"
b. Dead Soft 0.0016"
3. Subgingival extensions on select sizes for deeper preps
4. Full-Curve design to wrap further around tooth
5. Universal placement/removal tab

Firm Bands for superior strength and flexibility

Firm Bands are crafted from a proprietary stainless-steel alloy, offering exceptional strength and flexibility.

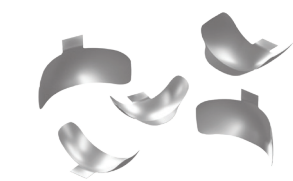

- 0.0015" thin!
- Resists deformation
 - Perfect for
 - Back-to-back
 - Deeper restorations
 - Unbroken contacts
 - Extra wide preps

Dead-Soft Bands for enhanced adaptation

The full curve design of these bands allows for optimal adaptation to the tooth surface, ensuring complete coverage of the preparation.



- 0.0016" thin!
- Easily customizable
 - Perfect for
 - Single Class IIs
 - General restorations